

---

# Experiment 10 Vinegar Analysis

---

*Experiment 10 Vinegar Analysis*

Downloaded from  
[evoluir.abar.org.br](http://evoluir.abar.org.br) by  
Matteo & Mathias

---

## INTRODUCTION TO EXPERIMENT 10 VINEGAR ANALYSIS

---

Vinegar is a staple in kitchens around the world, valued for its sharp flavor and preservative qualities. But beyond its culinary uses, vinegar presents a fascinating subject for chemical analysis.

**Experiment 10 Vinegar Analysis** is a classic laboratory procedure that introduces students and researchers to the principles of acid-base titration. This experiment is designed to determine the concentration of acetic acid in commercial vinegar, a measurement that has both practical and educational significance. Whether you are a chemistry student encountering titration for the first time or a curious enthusiast looking to understand the science behind everyday products, this comprehensive guide will walk you through every aspect of the analysis. By the end of this article, you will have a clear understanding of how to perform the experiment, interpret the results, and appreciate its real-world applications.

---

## THE PURPOSE OF EXPERIMENT 10 VINEGAR ANALYSIS

---

The primary goal of **Experiment 10 Vinegar Analysis** is to determine the percent by mass of acetic acid ( $\text{CH}_3\text{COOH}$ ) in a sample of vinegar. Acetic acid is the key component that gives vinegar its characteristic sour taste


and preservative properties. Commercial vinegar typically contains between 4% and 8% acetic acid by volume, though this can vary by type and brand. Through this experiment, students learn to apply titration techniques to quantify an unknown concentration, reinforcing concepts of molarity, stoichiometry, and solution chemistry. Moreover, the experiment cultivates essential lab skills such as precise measurement, careful observation, and accurate calculation.

---

## THEORETICAL BACKGROUND: ACID-BASE TITRATION

---

Titration is a quantitative analytical technique used to determine the concentration of an unknown solution by reacting it with a solution of known concentration. In **Experiment 10 Vinegar Analysis**, the unknown is vinegar (containing acetic acid), and the known solution (the titrant) is a standardized sodium hydroxide ( $\text{NaOH}$ ) solution. The reaction between acetic acid and sodium hydroxide is a classic neutralization reaction:

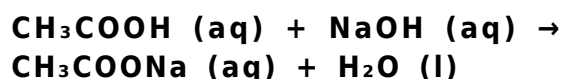

This reaction proceeds in a 1:1 molar ratio, meaning one mole of acetic acid reacts with one mole of sodium hydroxide. To detect the endpoint of the titration—the point at which the acid has been completely neutralized—a pH indicator is used. Phenolphthalein is the most common choice for this experiment because it changes from colorless to a faint pink color at the equivalence point

in a strong base-weak acid titration. The endpoint is visually distinct, allowing for accurate determination of the volume of titrant required.

## Why Standardization Matters

Before performing the vinegar analysis, the sodium hydroxide solution must be standardized. **Experiment 10 Vinegar Analysis** typically includes a preliminary step where NaOH is titrated against a primary standard, such as potassium hydrogen phthalate (KHP). This ensures that the concentration of the NaOH solution is precisely known, which is critical for accurate calculation of the acetic acid concentration. Without proper standardization, any error in the titrant concentration would propagate through the entire analysis, compromising the results.

### MATERIALS AND EQUIPMENT REQUIRED

To successfully conduct **Experiment 10 Vinegar Analysis**, you will need the following materials and equipment:

- **Vinegar sample:** White distilled vinegar is typically used, but any commercial vinegar can be analyzed.
- **Standardized sodium hydroxide solution:** Approximately 0.1 M, with a precisely known concentration.
- **Phenolphthalein indicator:** A 1% solution in alcohol.
- **Burette:** 50 mL capacity, for delivering the titrant.

- **Erlenmeyer flasks:** 250 mL, for containing the vinegar sample.
- **Graduated cylinder or volumetric pipette:** For measuring the vinegar sample volume.
- **Analytical balance:** For mass measurements if a mass-based approach is used.
- **Deionized water:** For diluting the sample and rinsing equipment.
- **Magnetic stirrer or stirring rod:** To ensure thorough mixing during titration.

Proper preparation of equipment is essential for accurate results. All glassware should be clean and dry, and the burette should be rinsed with the titrant solution before filling.

### PROCEDURE FOR EXPERIMENT 10 VINEGAR ANALYSIS

The following step-by-step procedure outlines the standard approach for **Experiment 10 Vinegar Analysis**. Always follow your instructor's specific guidelines and observe all safety protocols.

## Step 1: Standardization of Sodium Hydroxide (if required)

1. Weigh approximately 0.5 g of dried KHP accurately into an Erlenmeyer flask.
2. Dissolve the KHP in about 50 mL of

deionized water.

3. Add 2–3 drops of phenolphthalein indicator.
4. Titrate with the NaOH solution until a faint pink color persists for 30 seconds.
5. Record the volume of NaOH used and calculate the exact molarity.

## Step 2: Preparation of Vinegar Sample

1. Obtain a commercial vinegar sample. Note the label information for comparison.
2. Using a volumetric pipette, transfer exactly 5.00 mL of vinegar into a clean 250 mL Erlenmeyer flask.
3. Add approximately 50 mL of deionized water to dilute the sample. This dilution helps to reduce the intensity of the color change and improve endpoint detection.
4. Add 2–3 drops of phenolphthalein indicator.

## Step 3: Titration

1. Fill the burette with the standardized NaOH solution and record the initial volume.
2. Place the Erlenmeyer flask under the burette and begin adding NaOH slowly while swirling the flask constantly.
3. Watch for the appearance of a faint pink color. When the color starts to form but disappears with swirling, you are approaching the

endpoint.

4. Continue adding NaOH drop by drop until the faint pink color persists for at least 30 seconds.
5. Record the final volume of NaOH in the burette.

## Step 4: Replicates and Data Collection

Repeat the titration at least three times to ensure consistency. The volume of NaOH used should agree within  $\pm 0.1$  mL for precise results. Record all data in a well-organized table, including initial burette reading, final burette reading, and the volume of titrant used.

---

### DATA ANALYSIS AND CALCULATIONS

---

After collecting the titration data, the next phase of **Experiment 10 Vinegar Analysis** is calculation. The following steps will guide you from raw data to the final percentage of acetic acid in the vinegar sample.

## Calculating Moles of NaOH Used

For each titration trial, calculate the moles of NaOH that reacted:

**Moles of NaOH = Molarity of NaOH  
(mol/L)  $\times$  Volume of NaOH (L)**

# Determining Moles of Acetic Acid

## Example Calculation

### Acid

Because the reaction is 1:1, the moles of acetic acid are equal to the moles of NaOH at the endpoint:

**Moles of  $\text{CH}_3\text{COOH}$  = Moles of NaOH**

### Calculating the Mass of Acetic Acid

Convert moles of acetic acid to grams using its molar mass (60.05 g/mol):

**Mass of  $\text{CH}_3\text{COOH}$  = Moles of  $\text{CH}_3\text{COOH}$   $\times$  60.05 g/mol**

### Finding the Percent by Mass

Finally, calculate the percent by mass of acetic acid in the original vinegar sample:

**%  $\text{CH}_3\text{COOH}$  = (Mass of  $\text{CH}_3\text{COOH}$  / Mass of Vinegar Sample)  $\times$  100%**

The mass of the vinegar sample can be determined either by weighing the pipetted volume directly or by using the density of vinegar (typically 1.01 g/mL for white vinegar). If you used a 5.00 mL sample, the mass is approximately 5.05 g.

Suppose you used 0.105 M NaOH and titrated a 5.00 mL vinegar sample. The average volume of NaOH used in three trials is 18.50 mL. Moles of NaOH =  $0.105 \text{ mol/L} \times 0.01850 \text{ L} = 0.0019425 \text{ mol}$ . This equals the moles of acetic acid. Mass of acetic acid =  $0.0019425 \text{ mol} \times 60.05 \text{ g/mol} = 0.1166 \text{ g}$ . If the vinegar mass is 5.05 g, then % =  $(0.1166 \text{ g} / 5.05 \text{ g}) \times 100\% = 2.31\%$ . This value would be compared to the label claim, which is typically around 5% for white vinegar. Note that your results may differ based on dilution factors and sample volume adjustments.

---

### SOURCES OF ERROR AND HOW TO MINIMIZE THEM

---

No laboratory experiment is entirely free from error, and **Experiment 10 Vinegar Analysis** is no exception. Recognizing potential sources of error is a valuable part of the learning process. Common errors include:

- **Improper endpoint detection:** Adding too much titrant can overshoot the endpoint, leading to a false high reading. Slow, careful addition near the endpoint is crucial.
- **Parallax error:** Reading the meniscus at eye level is essential for accurate volume measurement.
- **Contaminated equipment:** Residue in the burette or flask can affect the reaction. Always rinse and clean thoroughly.
- **Incomplete standardization:** If the NaOH solution is not accurately standardized, all

subsequent calculations will be off.

- **Temperature fluctuations:**

Volume measurements can change with temperature. Ideally, perform the titration at room temperature.

To minimize errors, use calibrated equipment, practice careful technique, and perform multiple trials. Statistical analysis of your results, such as calculating the mean and standard deviation, provides insight into the precision of your work.

---

### SAFETY PRECAUTIONS FOR EXPERIMENT 10 VINEGAR ANALYSIS

---

Safety is paramount in any chemistry laboratory. Although vinegar and dilute sodium hydroxide are relatively safe, you should always follow standard lab safety protocols:

- Wear safety goggles at all times to protect your eyes from chemical splashes.
- Wear a lab coat or apron to protect your clothing.
- Handle sodium hydroxide carefully; it is caustic and can cause skin irritation.
- Avoid inhaling vinegar fumes,

especially from concentrated samples.

- Keep all chemicals away from open flames (though no flammable solvents are used in this experiment).
- Dispose of all solutions properly according to your laboratory's waste management guidelines.

---

### REAL-WORLD APPLICATIONS OF VINEGAR ANALYSIS

---

The skills and concepts learned in **Experiment 10 Vinegar Analysis** extend far beyond the classroom. The ability to determine the concentration of an acid in a consumer product has numerous practical applications:

- **Quality control in food manufacturing:** Companies must ensure that vinegar products meet regulatory standards for acetic acid content. The titration method used in this experiment is a standard quality control technique.
- **Home brewing and pickling:** Hobbyists and artisans often need to test acidity levels for safe preservation and consistent flavor.
- **Environmental monitoring:** Acid-base titration is used to measure the acidity of rainwater